


ROSA Ochratoxin Quantitative Test for Wine and Juice

**10
minutes**

Product Overview

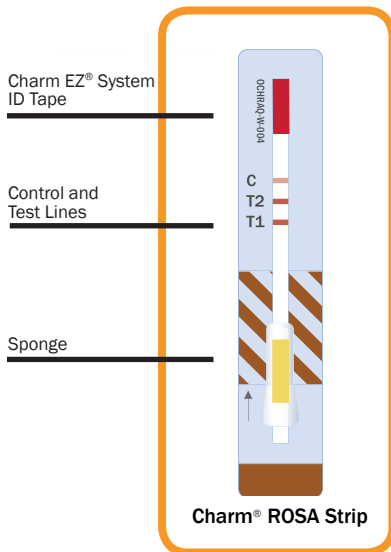

The **ROSA® Ochratoxin Quantitative Test for Wine and Juice (OCHRAQ-W)** delivers accurate results using ROSA (Rapid One Step Assay) lateral Flow technology. The diluted sample is applied to the test strip, incubated for 10 minutes and read in the Charm EZ-M system.


OHCRAQ-W Test:

- Range: 0 to 5000 ppt
- Results in 10 minutes

Run on the Charm EZ-M System:

- High volume output with batching using multiple quad incubators
- Standard curve built into reader reducing testing cost by eliminating end-user calibration
- Printable or downloadable results to any computer or LIMS system

Simple Procedure

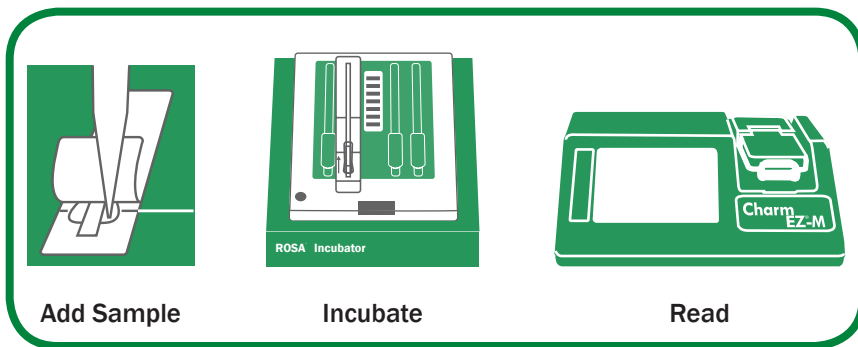

Charm Validated Commodities

Grape Juice	Wine		
-------------	------	--	--



659 Andover Street, Lawrence, MA 01843, USA
T +1.978.687.9200 | **F** +1.978.687.9216
E info@charm.com | **www.charm.com**

© 2016 Charm Sciences, Inc. Charm, Charm EZ and ROSA are registered trademarks, and test strip product packaging is trade dress, of Charm Sciences, Inc.

MRK-760-002 Feb-2016

Ordering Info

Order Codes

LF-OCHRAQ-W-20K
LF-OCHRAQ-W-100K
LF-OCHRAQ-W-500K

Each Kit Includes

Test Strips, Positive Control, Dilution Buffer

See Operator's Manual for additional kits and order codes