

# Field Incubator



## Charm Field Incubator for Portable Testing

The **Charm Field Incubator** is a portable incubator, making testing with Charm EZ Compatible<sup>®</sup> Strips fast and easy.

### Equipment Overview

The Charm<sup>®</sup> Field Incubator allows for easy compact incubation for all ROSA<sup>®</sup> (Rapid One Step Assay) tests. Powered by a rechargeable battery, the portable Field Incubator has software selection capabilities for temperature and incubation time. This self-contained unit can be operated on a truck, on the farm, or in the field.


### Features & Benefits

- ✓ Battery Operated – Lithium ion battery with 10-24 hour runtime
- ✓ Portable – Hand-held unit for use in lab or vehicle
- ✓ Rechargeable – Unit can be charged through wall adaptor
- ✓ Indicator Lights – Green and red indicators for instant visual of time and temperature
- ✓ Versatile – Can be used to run any Charm ROSA test
- ✓ Adjustable – Easy temperature and calibration adjustment


### Simple Procedure



Add Sample

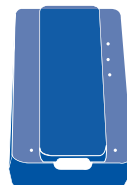

Incubate

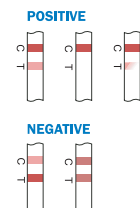


Visually Read

# Instrument Overview



**Ergonomic Design**

A durable molded casing allows for single-hand operation

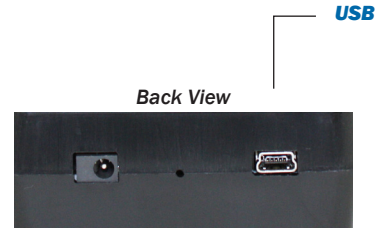

**Rechargeable / Replaceable Battery**  
24 hour runtime

**Portable Size**

**Software Configurable**

A pre-calibrated unit meets temperature and timer requirements for all Charm strips

# Specifications & Order Codes

Specifications

<b>Unit Dimensions</b>	7.3 x 4.2 x 11.5 cm (W x H X D)
<b>Weight</b>	290 g
<b>User Interface</b>	USB
<b>Power Source (AC)</b>	Lithium Ion Battery or 5 VDC
<b>Full Charge Battery Run Time</b>	Up to 24 hours
<b>Calibration (ATP)</b>	Factory Set (adjustments can be made through USB interface)

Order Info

Order Codes:	Contains:
LF-INC1-FIELD	Charm Field Incubator



659 Andover Street, Lawrence, MA 01843-1032, USA  
**T** +1.978.687.9200 | **F** +1.978.687.9216 | **E** info@charm.com | **www.charm.com**

© 2015 Charm Sciences, Inc. Charm, EZ Compatible, and ROSA are registered trademarks of Charm Sciences, Inc. Windows and Excel are registered trademarks of Microsoft Corporation in the United States and other countries.